

# Table of Contents

Nonlinear Curve Fitting: Fit Plot

1

*Further reading*

1

*Fit Function is a Sum of Fit Curves*

2

*Copying Fit Curves from One Fit Plot to Another*

2

*Setting Initial Values of Parameters*

2

*Parameter Locking*

2

*Fit Intervals*

2

*Baseline Fitting and Extraction*

3

*'Data-Baseline' Table Column*

3

*Viewing the Difference between Data and Fit Sum Function (Residual)*

4

*Fit One Curve*

4

*Joining the Parameters of Fit Curves*

4

*Fitting Algorithm*

4

*Fit Progress Window*

5

*Weighting of y data*

6

*Undoing Fit*

6

*Formulas*

6

*See Also*

7




# Nonlinear Curve Fitting: Fit Plot

Nonlinear least squares data fitting can be performed using Fit Plot. To create a Fit Plot, select x and y columns in table, then select **Table → Create Fit Plot** in main menu, or use context menu in table, or use Create Fit Plot button in the toolbar.

'Nonlinear' means here that analytical fitting function depends nonlinearly on varying parameters (fit parameters). Linear fitting is a quite simple method, which is based on solving the system of linear equations. Unlike linear fitting, nonlinear fitting is performed by iterative algorithm which needs the user to set the initial values of fit parameters.

To fit the data, implement these steps:

1. Create a Fit Plot, specify y weighting in Plot properties, if any
2. Specify fit function by adding Fit Curves
3. Specify initial values of fit parameters (drag curves or enter accurate values)
4. Specify used x data interval
5. Run fitting

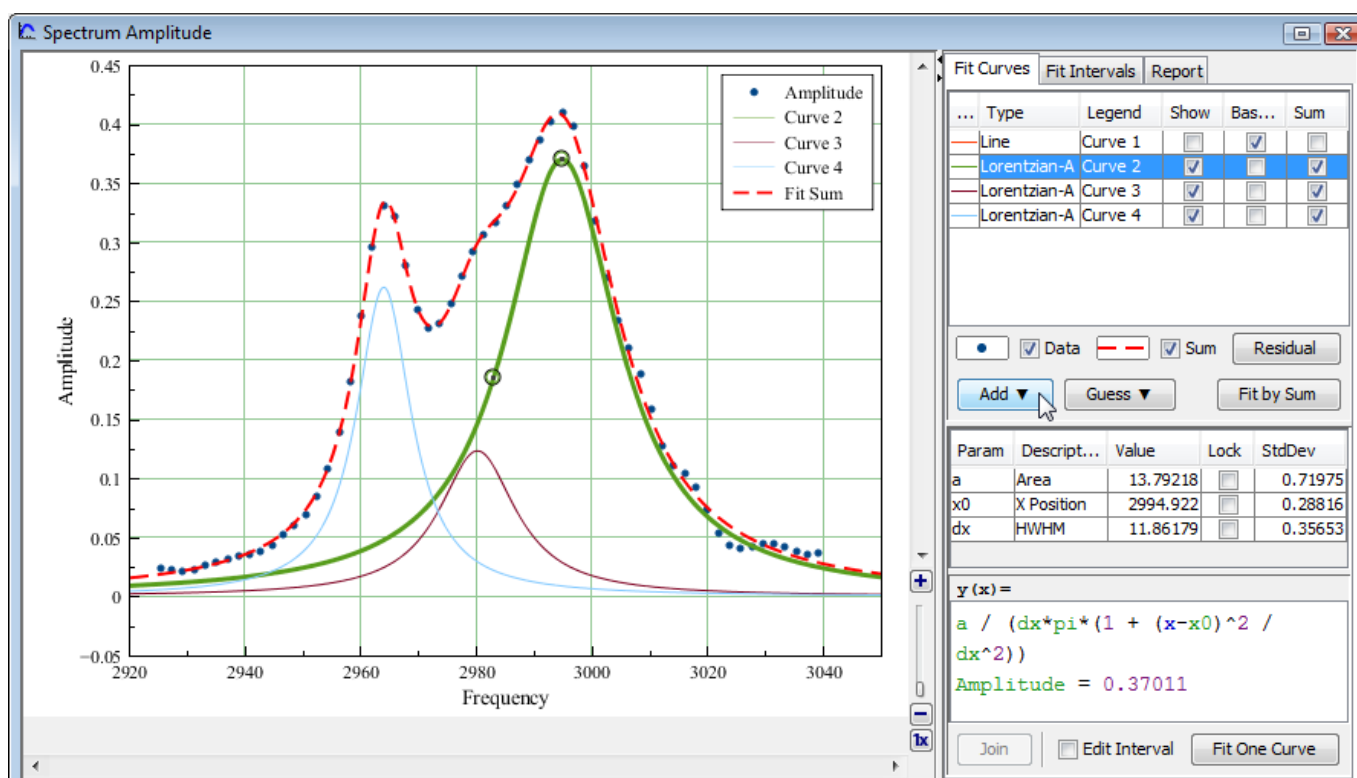

## Further reading

This manual does not completely cover the complex nonlinear fitting methodology. We recommend you to take a look at this book:

- H. Motulsky and A. Christopoulos, *Fitting Models to Biological Data Using Linear and Nonlinear Regression: A Practical Guide to Curve Fitting*. 2003, GraphPad Software Inc., San Diego CA, graphpad.com. PDF is available [here](#).

## Fit Function is a Sum of Fit Curves

MagicPlot considers fit function as a **sum** of Fit Curves. Ordinarily in peaks fitting each Fit Curve corresponds to one peak in experimental data. There is a number of predefined Fit Curves (Line, Parabola, Gauss, Lorentz, etc.) You can also specify a custom Fit Curve. Baseline fitting components may be added to the fitting sum, too.

Fit Plot window contains the list of Fit Curves. Each Fit Curve in the list has three checkboxes: Show, Baseline, Sum:

- **Show:** Specifies whether to show this Fit Curve on plot. Active only if Baseline checkbox is not set
- **Baseline:** Toggles the subtracting of this Fit Curve from experimental data. You also can use **Residual** button to subtract Fit Sum from data
- **Sum:** Specifies whether to use this Fit Curve in sum fit function

Below the Fit Curves list is a parameters table which shows names, values, and descriptions of parameters relating to selected Fit Curve.

## Copying Fit Curves from One Fit Plot to Another

You can copy and paste Fit Curves in curves table as usual. Use context menu in curves table or press Ctrl+C, Ctrl+V on PC and Cmd C, Cmd V on Mac (curves table must have focus, not the plot itself). MagicPlot copies Fit Curves parameters values, baseline/sum/show checkbox values and style.

## Setting Initial Values of Parameters

Nonlinear fitting assumes that certain initial values of parameters are set before fitting. This procedure is very easy if you use predefined Fit Curves: you can drag curves on plot.

Initial parameters values for each Fit Curve can also be set in parameter table.

## Parameter Locking


You can lock parameter(s) to prevent varying this parameter during fit and to prevent its changing due to setting initial values by mouse dragging (for built-in functions). Set the checkbox in Lock column in parameters list.

## Fit Intervals

You can set the x intervals of the data. Data points outside these intervals are not used to compute the minimizing residual sum of squares (see below). You can use this feature if some data points (especially in the beginning or the end) are inaccurate, e.g. noisy.

Select **Fit Interval** tab to set intervals visually or edit accurate borders values:

- Double click on interval to split it
- Drag the interval border to move it. If intervals intersect, they will be merged
- Use context menu on the plot to create, delete and split intervals

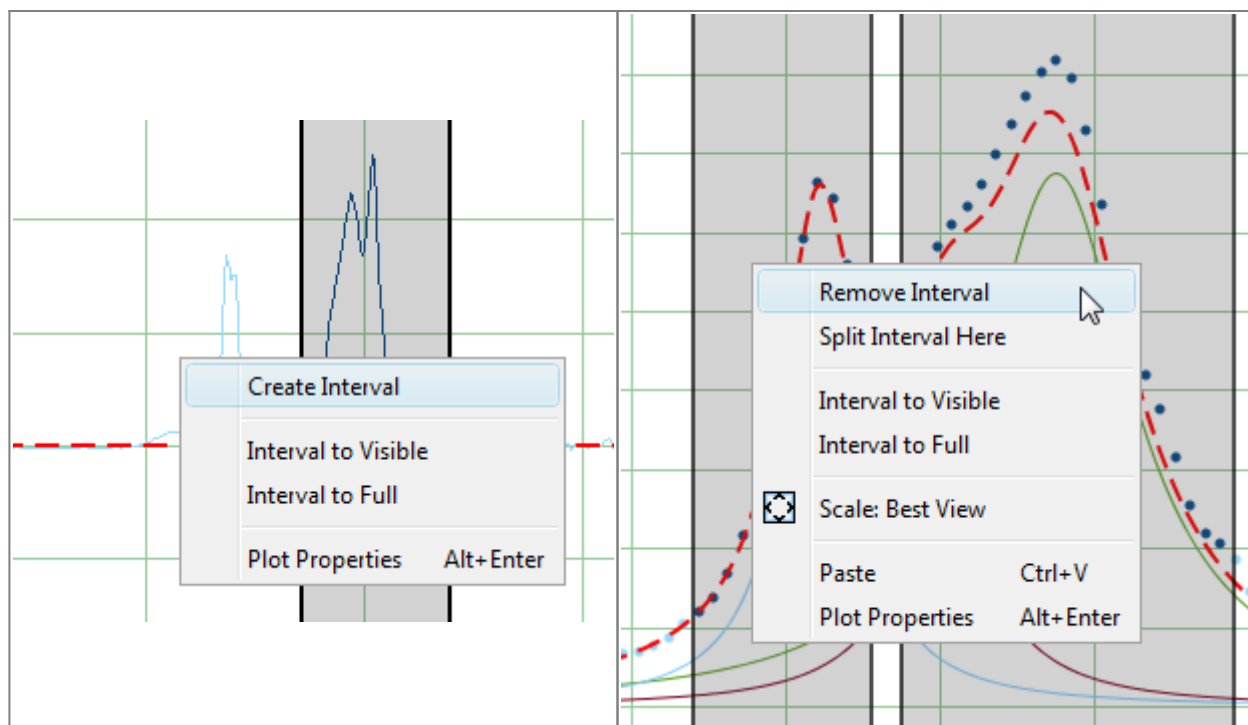

## Baseline Fitting and Extraction

Fit Interval is also usable when baseline fitting. Before baseline fitting you can specify the interval which does not contain any signal points and contains baseline only. Set **Baseline** checkboxes at baseline Fit Curves after baseline fitting to subtract baseline from data. Then specify the whole interval and fit the data.

The most appropriate curve type for baseline fitting is [spline](#).

Note that if you execute one of data processing algorithms (integration, FFT, etc.) on Fit Plot, then the difference between the data and baseline curves (which you do see on the plot) will be processed. You can use this behaviour to exclude baseline from data before integrating, see [Integration](#) for more information.

## 'Data-Baseline' Table Column

The 'Data-Baseline' column is appended to the Table with initial (x, y) fit data when you create Fit Plot. The 'Data-Baseline' column contains the difference between initial y data and baseline approximation (the sum of Fit Curves for which **Baseline** checkbox is set).

It is 'Data-Baseline' column that is actually plotted on Fit Plot.

Use 'Data-Baseline' column in Table if you want to process the data without baseline. This column is also used as initial data if you use **Processing** menu when Fit Plot is active.

## Viewing the Difference between Data and Fit Sum Function (Residual)

MagicPlot offers two different ways to view the difference between data and Fit Sum function:

- You can set Baseline checkboxes for all fit function components to subtract them from data and explore the residual plot

or

- You can press and hold Residual button. The difference between data and Fit Sum function is shown while button is pressed. You can use either mouse or space key (if button is selected) to hold Residual button.

## Fit One Curve

You can also use MagicPlot to fit the data with single selected Fit Curve by pressing Fit One Curve button. In this case a specific data interval for each Fit Curve is used and the main fitting data interval (set in Fit Interval tab) is ignored. Select Set Interval checkbox in the bottom of the Fit Plot panel to set specific fit intervals for each Fit Curve.

Because of using individual data interval this method is useful for baseline fitting. In order to fit baseline specify the intervals which does not contain signal (peaks) and contain only noise.

## Joining the Parameters of Fit Curves

MagicPlot allows coupling of fit parameters. See [Joining the Parameters of Fit Curves](#) for details.

## Fitting Algorithm

MagicPlot uses iterative [Levenberg-Marquardt](#) [nonlinear least squares](#) curve fitting algorithm which is widely used in most software.

Fit procedure iteratively varies the parameters  $\beta_k$  of fit function  $f(x, \beta_1, \dots, \beta_p)$  to minimize the residual sum of squares (RSS,  $\chi^2$ ):

here:

$$\chi^2 = \sum_{i=1}^N \left( w_i (y_i - f(x_i, \beta_1, \dots, \beta_p)) \right)^2 \rightarrow \min$$

- $x_i$  and  $y_i$  are the data points,
- $N$  is total number of points,
- $f(x, \beta_1, \dots, \beta_p)$  is the fit function which depends on value of  $x$  and fit parameters  $\beta_k$ ,
- $p$  is the number of fit parameters  $\beta_k$ ,
- $w_i$  are normalized  $y$  data weighting coefficients for each point  $y_i$ :

Calculation of the new guess of parameters on each fit iteration is based on the fit function partial derivatives for current values of fit parameters and for each  $x$  value:

$$\frac{\partial f}{\partial \beta_m} \left( x_i, \beta_1, \dots, \beta_p \right)$$

To start minimization, you have to provide an initial guess for the parameters.

## Fit Procedure Stop Criteria

After each iteration except the first MagicPlot evaluates *deviation decrement*  $D$ :

$$D = \frac{\sum_{i=1}^N x_i^2}{\sum_{i=1}^N x_i^2}$$

Deviation decrement shows how the residual sum of squares (RSS) on current iteration relatively differs from that on the previous iteration.

The iterative fit procedure stops on one of two conditions:

- If the deviation decrement  $D$  is less than minimum allowable deviation decrement, which is  $10^{-9}$  by default


or (and)

- If the number of iterations exceeds maximum number of iterations, which is 100 by default

You can change the minimum allowable deviation decrement and maximum number of iterations in Fitting tab of MagicPlot Preferences.

## Fit Progress Window

MagicPlot indicates fit process with a special window. Fitting curves are periodically updated on plot while fitting so you can see how fit converges.

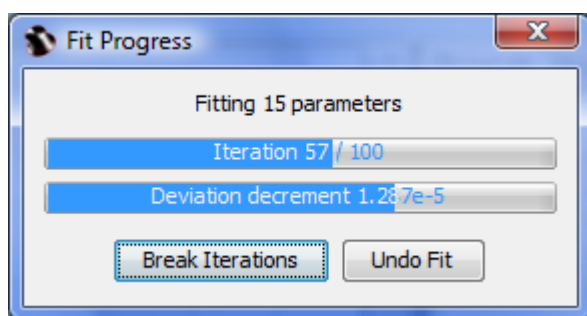

MagicPlot shows current iteration number and deviation decrement with two progress bars while fit is performed. The fit process stops when one of these progress bars reaches the end.

You can see two buttons on fit progress window:

- **Break Iterations:** Breaks iterations after current iteration. Use this button if you suspect that further iterations will not change the result.
- **Undo Fit:** Breaks iterations and reverts fit parameters to their initial (before fit) values. Use this button if you see that fit process converges to wrong result; change initial values of parameters and run fit again.

# Weighting of y data

MagicPlot can use weighting of y values based on y errors  $s_i$ :

- If standard y errors are **not** specified: all  $w_i=1$
- If standard y errors  $s_i$  are specified:

$w_i = C \frac{1}{s_i^2}$   
here C is normalizing coefficient (to make the sum of  $w_i$  be equal to one):

$C = \frac{1}{\sum_{i=1}^N s_i^2}$   
In Fit Plot Properties dialog (Plot Data tab) you can set one of the following methods to evaluate standard y errors  $s_i$ :

- Get y errors from table column(s),
- Percent of data for every point,
- Fixed value or Standard deviation — do not use in weighting because in this case the error values are the same for all data points.

# Undoing Fit

You can undo fit and undo changing initial parameters as usual using Undo function. It is a handy feature when experimenting with different models and initial parameters (e.g. peaks positions).

# Formulas

In the table below you can find the formulas which MagicPlot uses to calculate fit parameters and values in Fit Report tab.

Because of some confusion in the names of the parameters in different sources (books and software), we also give many different names of same parameter in *note* column.

| Parameter Name                            | Symbol                    | Formula                                      | Note  |
|---|---------------------------|--|---|
| Original Data and Fit Model Properties    |                           |  |   |
| Number of used data points                | $N$                       | —  | This is the number of data points inside specified Fit Interval.  |
| Fit parameters                            | $\beta_1, \dots, \beta_p$ | —  | For peak-like functions (Gauss, Lorentz) these parameters are amplitude, position and half width at half maximum. |
| Number of fit function parameters $\beta$ | $p$                       | —  | This is the total number of parameters of all fit curves which are summarized to fit.                             |
| Degrees of freedom                        | $\nu$                     | $N - p$                                      |   |
| Estimated mean of data                    | $\bar{y}$                 | $\frac{1}{N} \sum_{i=1}^N y_i$               |   |
| Estimated variance of data                | $s^2$                     | $\frac{1}{N} \sum_{i=1}^N (y_i - \bar{y})^2$ | Not used by fit algorithm, only for comparison.   |
| Data total sum of squares, TSS            | TSS                       | $\sum_{i=1}^N (y_i - \bar{y})^2$             | TSS is also called sum of squares about the mean and acronym SST is also used.                                    |



| Parameter Name   | Symbol           | Formula   | Note   |
|--|------------------|---|--|
| Original Data and Fit Model Properties                 |                  |   |  |
| Fit Result   |                  |   |  |
| Residual sum of squares, RSS                           | $\chi^2$         | $\sum_{i=1}^N \left( y_i - f(x_i, \beta_1, \dots, \beta_p) \right)^2$ | This value is minimized during the fit to find the optimal fit function parameters.<br><i>RSS is also called the sum of squared residuals (SSR), the error sum of squares (ESS), the sum of squares due to error (SSE).</i>  |
| Reduced $\chi^2$                                       | $\chi^2_{red.}$  | $\frac{\chi^2}{N-p}$  | The advantage of the reduced chi-squared is that it already normalizes for the number of data points and model (fit function) complexity.<br><i>Reduced <math>\chi^2</math> is also called mean square error (MSE) or the residual mean square.</i>  |
| Residual Standard Deviation                            | $s$              | $\sqrt{\frac{\chi^2_{red.}}{N-p}}$                                    | Standard deviation is also called root mean square of the error (Root MSE)   |
| Coefficient of determination                           | $R^2$            | $\frac{\chi^2_{red.}}{1 - \frac{\chi^2_{red.}}{TSS}}$                 | $R^2$ will be equal to one if fit is perfect, and to zero otherwise. This is a biased estimate of the population $R^2$ , and will never decrease if additional fit parameters (fit curves) are added, even if they are irrelevant.   |
| Adjusted $R^2$   | $\overline{R}^2$ | $1 - \frac{N-1}{N-p-1} (1 - R^2)$                                     | Adjusted $R^2$ (or <i>degrees of freedom adjusted R-square</i> ) is a slightly modified version of $R^2$ , designed to penalize for the excess number of fit parameters (fit curves) which do not add to the explanatory power of the regression. This statistic is always smaller than $R^2$ , can decrease as you add new fit curves, and even be negative for poorly fitting models |
| Covariance Matrix of Parameters $\beta_k$              | $\Sigma$         | $\frac{\chi^2}{N-p} \alpha^{-1}$                                      | Here $\alpha$ is the matrix of partial derivatives of fit function with respect to parameters $\beta_m$ and $\beta_n$ which is used for fitting:   |
| Standard Deviation of Parameters $\beta_k$ , std. dev. | $s_k$            | $\sqrt{\Sigma_{k,k}}$   | $\alpha_{m,n} = \sum_{i=1}^N w_i \frac{\partial f}{\partial \beta_m} \left( x_i, \beta_1, \dots, \beta_p \right) \frac{\partial f}{\partial \beta_n} \left( x_i, \beta_1, \dots, \beta_p \right)$  |
| Correlation Matrix of Parameters $\beta_k$             | $C$              | $C_{i,j} = \frac{\Sigma_{i,j}}{s_i s_j}$                              |  |

See Also

- [Specifying Custom Fit Equation](#)
- [Using Spline for Baseline Subtraction](#)
- [Joining the Parameters of Fit Curves](#)
- [Guessing Peaks](#)
- [Predefined Fit Curves Equations](#)
- [transform\\_xy](#)
- [Calculating Integrals and Statistics on Intervals using Fit Plot](#)

From:  
<https://magicplot.com/wiki/> - **MagicPlot Manual**

Permanent link:  
<https://magicplot.com/wiki/fitting?rev=1298370371>

Last update: **Sun Nov 8 12:20:32 2015**

